

## Platinum Prosthesis for Stapedectomy

CODE  
SPL 03.00

This prosthesis is made entirely of 999‰ fine platinum ribbon.


Platinum is a highly malleable and well tolerated metal.

The hook has three transverse grooves: under slight pressure of a small forceps it bends only in the three points corresponding to the grooves

gently wrapping the long process of the incus without causing ischemia and without the danger of anomalous closures.

The lower extremity of the stem is elbow shaped to provide anchorage to the new membrane of the oval window.

The weight of the prosthesis (approximately 3 mg., with variations depending on the size) is very close to the weight of the stapes.

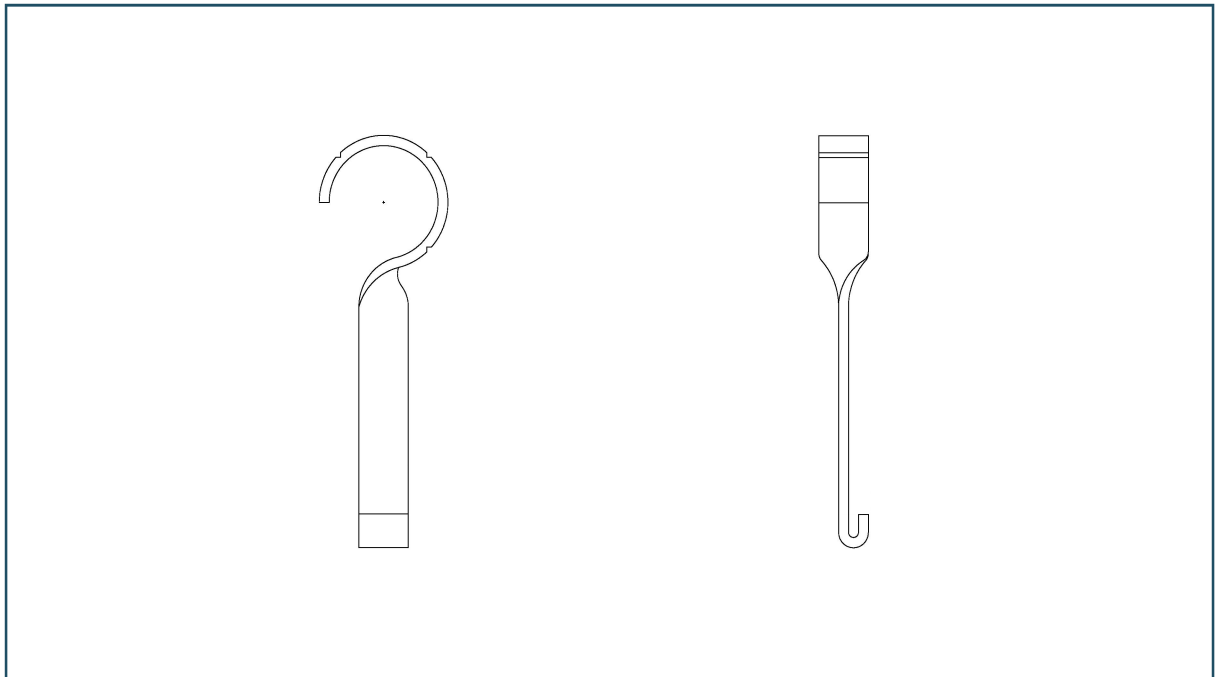


Design: C.Zini, Parma, Italy (using Scola's original design)

**Dimensions**

Total length (mm):  
3.00 - 3.50 - 4.00  
4.25 - 4.50 - 4.75  
5.00 - 5.50 - 6.00

**Packaging**

Dispenser box containing 6  
sterile, individually packaged  
implants.

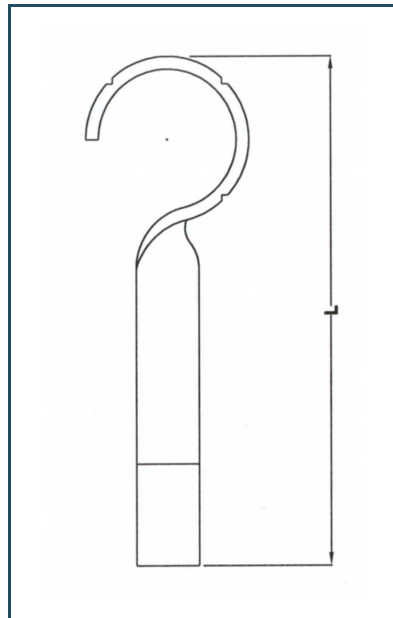

Total length (L) mm	Code
3.00	SPL 03.00.300
3.50	SPL 03.00.350
4.00	SPL 03.00.400
4.25	SPL 03.00.425
4.50	SPL 03.00.450
4.75	SPL 03.00.475
5.00	SPL 03.00.500
5.50	SPL 03.00.550
6.00	SPL 03.00.600

