

Platinum and Audio-PE Prosthesis for Stapedectomy

CODE
SPL 03.04

PE: High molecular weight porous polyethylene with pores of approximately 100 microns
This material stimulates the ingrowth of connective tissue and blood vessels


This prosthesis is made entirely of 999‰
fine platinum ribbon.

Platinum is a highly malleable and well
tolerated metal.

The hook has three transverse grooves:
under slight pressure of a small forceps
it bends only in the three points corresponding
to the grooves

gently wrapping the long process of the incus
without causing ischemia
and without the danger of anomalous
closures.

The AUDIO-PE® distal end provides excellent
anchorage to the new membrane of the oval
window due to its microporous structure
which encourages the ingrowth of
connective tissue and blood vessels.

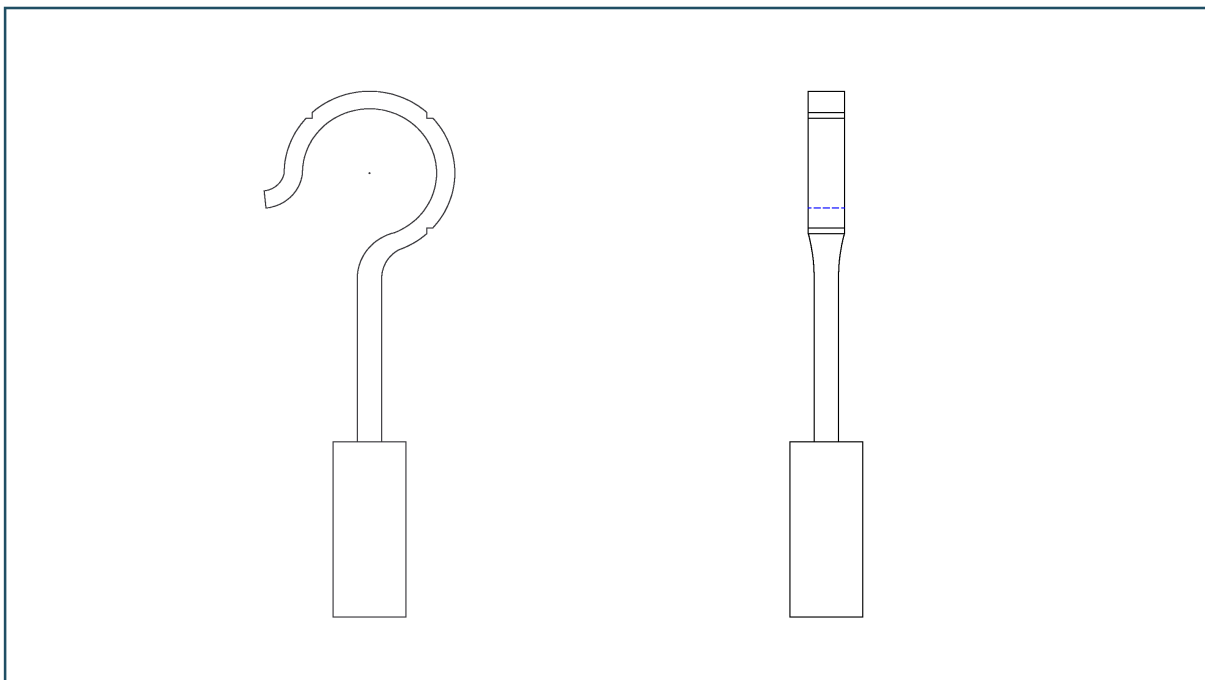

Dimensions

Total length (mm):

3.00 - 3.50 - 4.00


4.25 - 4.50 - 4.75

5.00 - 5.50 - 6.00

Piston diameter: 0.60 mm.

Packaging

Dispenser box containing 6 sterile, individually packaged implants.

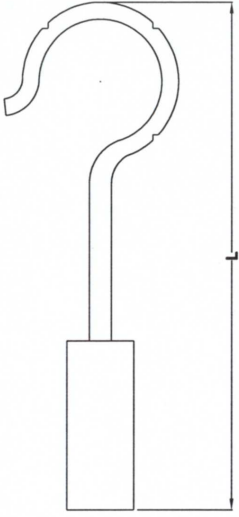	Total length (L) mm.	Code
	3.00	SPL 03.04.300
	3.50	SPL 03.04.350
	4.00	SPL 03.04.400
	4.25	SPL 03.04.425
	4.50	SPL 03.04.450
	4.75	SPL 03.04.475
	5.00	SPL 03.04.500
	5.50	SPL 03.04.550
6.00	SPL 03.04.600	

